



HPPS/HT Filters

High Temperature Dual Layered
PES Membrane



For sterile applications at elevated temperatures (up to 203°F), HPPS/HT filters remove high particle loads from aqueous liquids. These validated dual layered PES cartridge filters are available in pore sizes that range from 0.03 μm to 1.0 μm .


Designed with a High Capacity PES prefilter layer and PES final filter layer, plus high temperature polypropylene membrane support and cartridge hardware, HPPS/HT filters deliver high flow across a wide pH range. The HPPS/HT filters' low binding characteristics make them ideal for filtering products with preservatives and proteins that can adsorb to media. They are flushed to remove manufacturing debris and reduce extractables. Products are 100% integrity tested.

Critical Process provides unrivaled delivery times, technical consulting before purchasing, and very competitively priced high-performance products. Our comprehensive testing & analysis and validation services support your team whenever they need it. Your process experts partnering with our filtration experts is how we deliver your company's solution right the first time.

Sterilizing Filters



CARTRIDGES – Nominal Dimensions
Length: 5 to 40 in. (12.7 to 101.6 cm)
Outside Diameter: 2.75 in. (7.0 cm)

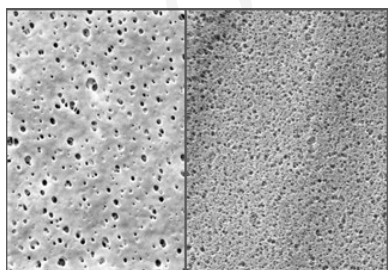

High Capacity PES / Asymmetric PES

HPPS/HT is recommended for:

- SVPs & LVPs
- Diagnostics
- Buffers
- WFI, Water Purification
- Vaccines
- Biologicals
- Ophthalmics

Maximum Operating Parameters

CARTRIDGES	
Operating Temperature (water)	203 °F at 30 psid (95 °C at 2.07 bard)
Forward Differential Pressure	80 psid at 68 °F (5.52 bard at 20 °C)
Reverse Differential Pressure	50 psid at 68 °F (3.45 bard at 20 °C)
Recommended Changeout Pressure	35 psid (2.41 bard)

Sanitization & Sterilization

Filtered Hot Water*	90 °C (194 °F), 30 minutes, multiple cycles, max 3 psid forward flow
Inline Steam*	275 °F (135 °C), 30 min, 25+ cycles
Autoclave*	250 °F (121 °C), 30 min, 25+ cycles
Chemical Sanitization	Performed using industry standard concentrations of hydrogen peroxide, peracetic acid, sodium hypochlorite and other selected chemicals.

*For all elevated temperature procedures above, a stainless-steel support ring is required.

Filtration Area (Nominal)

CARTRIDGES					
Length	5"	10"	20"	30"	40"
	12.7cm	25.4cm	50.8cm	76.2cm	101.6cm
Area	2.5 ft ²	5.4 ft ²	10.8 ft ²	16.2 ft ²	21.6 ft ²
	0.23m ²	0.50m ²	1.00m ²	1.51m ²	2.01m ²

Integrity Testing

PORE SIZE			DIFFUSION TEST PRESSURE						
μm	PSIG	BARG	DIFFUSION SPECIFICATIONS* (Final Layer Pore Size)						
			Length	2"	5"	10"	20"	30"	40"
0.1/0.03	60	4.14	mL/min (0.03μm, 0.10μm)	≤ 2.9	≤ 8.4	≤ 20	≤ 40	≤ 60	≤ 80
0.2/0.10	48	3.30							
0.5/0.22	35	2.41	mL/min (All Other Pore Sizes)	≤ 2.1	≤ 6.3	≤ 15	≤ 30	≤ 45	≤ 60
1.0/0.22	35	2.41							
1.0/0.45	20	1.37							
1.0/0.65	15	1.03							
1.0/0.80	12	0.82							
1.0/1.0	8	0.55							

* For water-wetted membrane

Construction Materials

Filtration Media	High Capacity PES membrane on polyester support prefilter layer and PES membrane final filter layer
Media Support	High Temperature Polypropylene
End Caps, Center Core, Outer Support Cage	High Temperature Polypropylene
Sealing Method	Thermal Bonding
O-Rings/Gaskets	Buna, Viton® (or FKM), EPDM, Silicone, FEP Encapsulated Silicone, FEP Encapsulated Viton (or FKM)

Validation

HPPS/HT filters are validated using test procedures that comply with ASTM F 838-15(ae1) protocols for the determination of bacterial retention in filters used for liquid filtration. The challenge level is a minimum of 10^7 organisms per cm^2 of filter media. CPF filters have > 7-log removal when challenged with the organisms listed below (0.03 μm , 0.10 μm and 0.22 μm meet the FDA definition of sterilizing grade filters).

0.03 μm : *Acholeplasma laidlawii*

0.10 μm : *Brevundimonas diminuta*

0.22 μm : *Brevundimonas diminuta*

0.45 μm : *Serratia marcescens*

0.65 μm : *Saccharomyces cerevisiae*

Endotoxins

The levels of bacterial endotoxins in aqueous extracts from HPPS/HT filters are below current USP limits as specified for water for injection.

Extractables

HPPS/HT filters typically exhibit low levels of non-volatile residues.

TOC and Conductivity

HPPS/HT filters conform with TOC standards of USP <643> and the water conductivity standards of USP <645> after an appropriate flush with purified water.

Toxicity Compliance

Materials used to construct HPPS/HT filters are non-toxic and meet the requirements for the MEM Elution Cytotoxicity Test and the requirements for Biological Reactivity Tests in the current version of the United States Pharmacopeia (USP) for Class VI - 121 °C Plastics.



Non-Fiber Releasing

HPPS/HT filters comply with Title 21 CFR sections 210.3 (b)(6) and 211.72, for non-fiber releasing filters.

FDA Compliance

Materials meet the requirements listed by the FDA as appropriate for use in articles intended for repeated food contact as specified in Title 21 CFR sections 174.5, 177.1500, 177.1520, 177.1630, 177.2440, and 177.2600 as applicable.

Flow Rates for HPPS/HT Cartridges by Pore Size

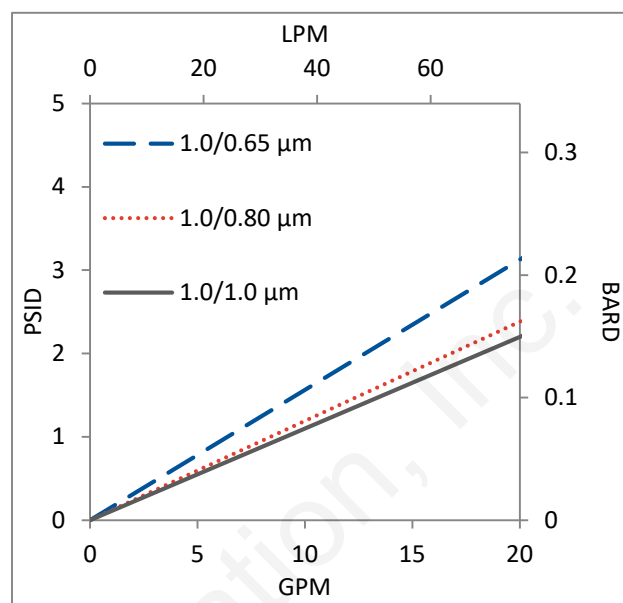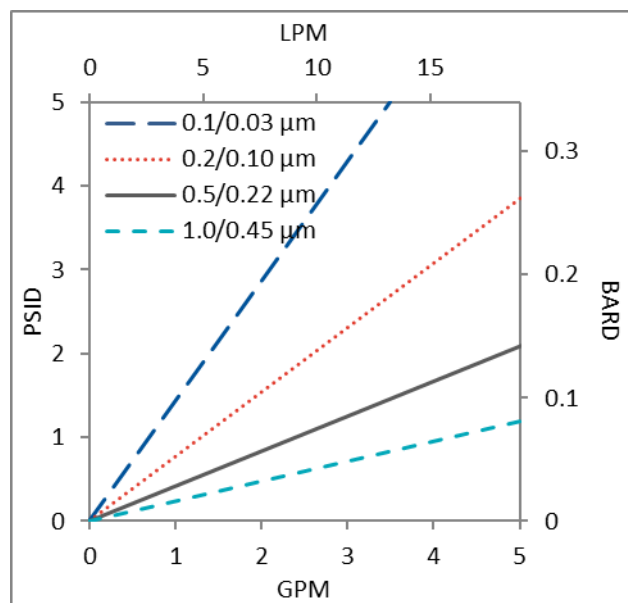

Flow rates for Cartridge filters are per 10-inch length. The test fluid is water at ambient temperature.


HPPS/HT Filters Ordering Information

Fill in the corresponding codes in the boxes below to build your Part Number.

To consult with one of our technical team members, request a quote or place an order:

call (603) 880-4420 or [contact us here](#).

Cartridge Filters

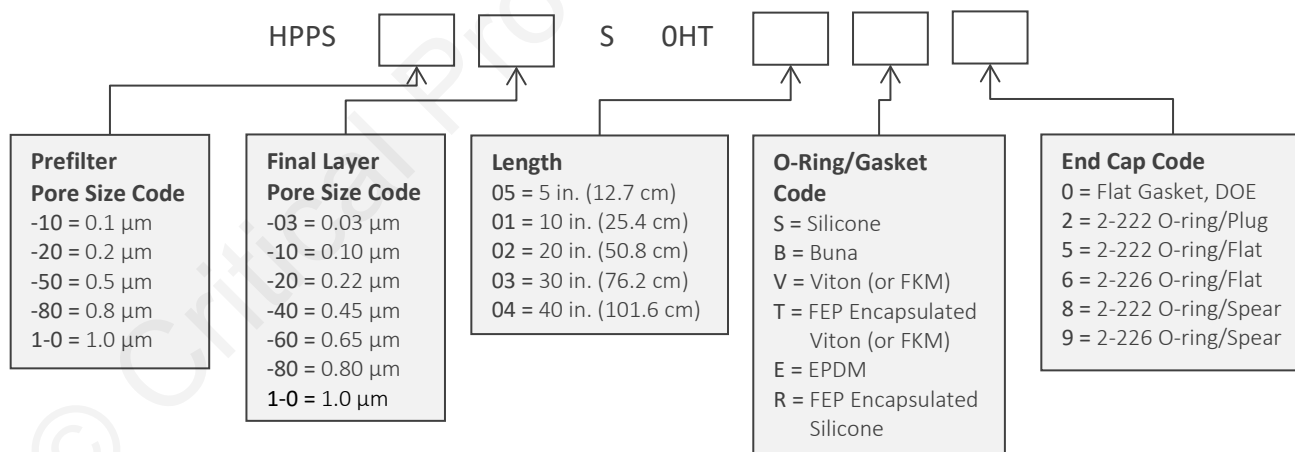

One Chestnut Street
Nashua, NH 03060
603.880.4420
FAX: 603.880.4536
CriticalProcess.com

The information contained herein is subject to change without notice. The Critical Process Filtration logo is a trademark of Critical Process Filtration, Inc. Viton is a trademark of DuPont Performance Elastomers L.L.C.
© 2021 Critical Process Filtration, Inc. • All Rights Reserved

Data Sheet HPPS/HTDS Rev C