



DPPC Micro Capsule Filters

Positively Charged, Dual Layered
PES Membrane

DPPC Micro Capsule filters consist of two layers of positively charged Polyethersulfone (PES) membrane.

The bioburden reduction prefilter and the sterilizing grade final filter each come in several pore sizes so you can configure to your system's needs. Available in 0.03, 0.10, 0.22, and 0.45 μm , DPPC Micro Capsules are validated for absolute bacteria retention to provide reliable sterile filtration performance. The positive charge removes negatively charged biological contaminants such as endotoxins, viruses, and other cell fragments. Depending on the level of contaminant and flow rate, DPPC filters will typically exhibit > 3-log removal of endotoxin.

This combination of functionality makes the DPPC filter an excellent choice for pharmaceutical and biopharmaceutical applications.

Critical Process provides unrivaled delivery times, technical consulting before purchasing, and competitively priced high-performance products. Our comprehensive testing & analysis, and validation services support your team whenever they need it. Your process experts partnering with our filtration experts is how we deliver your company's solution right the first time.

Sterilizing Filters

Endotoxin Removal




MICRO CAPSULES – Nominal Dimensions

Body Length: 1.9 in. (4.8 cm)

Overall Length: 2.8 to 3.8 in. (7.1 to 9.7 cm)

Outside Diameter: 2.6 in. (6.6 cm)

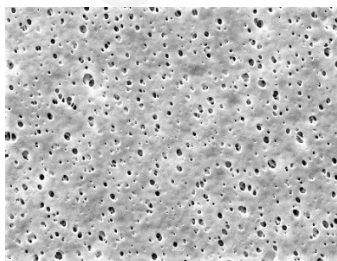

DPPC Micro Capsule filters are recommended for:

- Water Purification
- Water for Injection (WFI)
- Serum
- Plasma
- Vitamins

Maximum Operating Parameters

MICRO CAPSULES	
Liquid Operational Pressure	80 psi at 68 °F (5.52 bard at 20 °C)
Gases Operational Pressure	60 psi at 68 °F (4.14 bar at 20 °C)
Operating Temperature (water)	110 °F at 30 psid (43 °C at 2.07 bard)
Forward Differential Pressure	50 psid at 68 °F (3.45 bard at 20 °C)
Reverse Differential Pressure	40 psid at 68 °F (2.76 bard at 20 °C)
Recommended Changeout Pressure	35 psid (2.41 bard)

Sanitization & Sterilization

Autoclave	250 °F (121 °C), 30 min, 5+ cycles
Chemical Sanitization	Performed using industry standard concentrations of hydrogen peroxide, peracetic acid, sodium hypochlorite and other selected chemicals.

Integrity Testing

PORE SIZE	BUBBLE POINT MINIMUM*	
	PSIG	BARG
μm		
0.03	**	**
0.10	**	**
0.22	50	3.4
0.45	25	1.7

* For water wetted membrane

** Actual bubble point exceeds operational limits of Micro capsule filters.

Filtration Area (Nominal)

Area	0.52 ft ²
	483 cm ²

Construction Materials

Filtration Media	Positively Charged Double Layer PES Membrane
Media Support	Polypropylene
End Caps, Center Core, Outer Support Cage, Micro Capsule Housing	Polypropylene
Sealing Method	Thermal Bonding

Validation

DPPC Micro Capsule filters are validated using test procedures that comply with ASTM F 838-15(ae1) protocols for the determination of bacterial retention in filters used for liquid filtration. The challenge level is a minimum of 10^7 organisms per cm^2 of filter media. CPF filters have > 7-log removal when challenged with the organisms listed below (0.03 μm , 0.10 μm and 0.22 μm meet the FDA definition of sterilizing grade filters).

0.03 μm : *Acholeplasma laidlawii*

0.10 μm : *Brevundimonas diminuta*

0.22 μm : *Brevundimonas diminuta*

0.45 μm : *Serratia marcescens*

Endotoxins

The levels of bacterial endotoxins in aqueous extracts from DPPC Micro Capsule filters are below current USP limits as specified for water for injection.

Extractables

DPPC Micro Capsule filters typically exhibit low levels of non-volatile residues.

TOC and Conductivity

DPPC Micro Capsule filters conform with TOC standards of USP <643> and the water conductivity standards of USP <645> after an appropriate flush with purified water.

Toxicity Compliance

Materials used to construct DPPC Micro Capsule filters are non-toxic and meet the requirements for the MEM Elution Cytotoxicity Test and the requirements for Biological Reactivity Tests in the current version of the United States Pharmacopeia (USP) for Class VI - 121 °C Plastics.


Non-Fiber Releasing

DPPC Micro Capsule filters comply with Title 21 CFR sections 210.3 (b)(6) and 211.72, for non-fiber releasing filters.

FDA Compliance

Materials meet the requirements listed by the FDA as appropriate for use in articles intended for repeated food contact as specified in Title 21 CFR sections 174.5, 177.1500, 177.1520, 177.1630, 177.2440, and 177.2600 as applicable.

Flow Rates of Final Layer for DPPS Micro Capsule by Pore Size

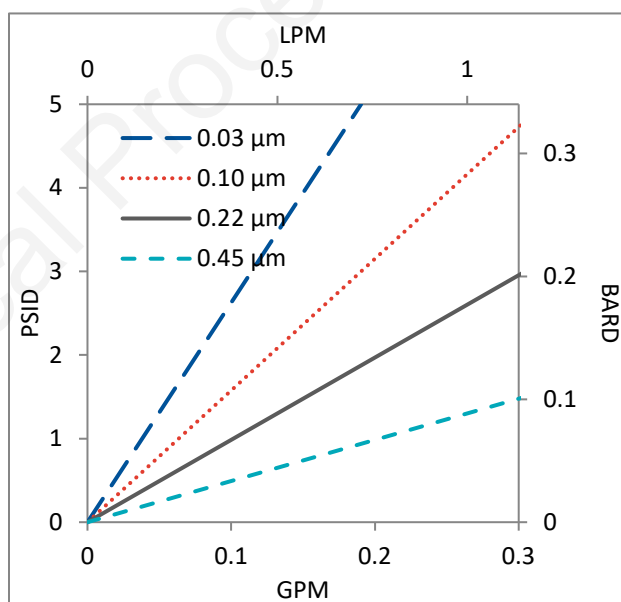

Flow rates for Micro Capsule filters are per filter. The test fluid is water at ambient temperature. Flows are tested using a Micro capsule filter with 1/2" Sanitary inlet and outlet ports. Rates will vary based on end configuration of the Micro capsule.

DPPC Micro Capsule Filters Ordering Information

All Critical Process filters are configurable to meet customer specifications.
Fill in the corresponding codes in the boxes below to build your Part Number.

To consult with one of our technical team members, request a quote or place an order:
call (603) 880-4420 or [contact us here](#).

****Please note this product is not designed or approved for use in Hemodialysis applications****

Micro Capsule Filters

MIC DPPC 000 -

Prefilter Pore Size Code	Final Layer Pore Size Code	Inlet	Inlet Vent Port	Outlet	Outlet Vent Port	Side Vent O-Ring**
-10 = 0.10 µm -20 = 0.22 µm -40 = 0.45 µm -60 = 0.65 µm -80 = 0.80 µm 1-0 = 1.0 µm 1-2 = 1.2 µm	-03 = 0.03 µm -10 = 0.10 µm -20 = 0.22 µm -40 = 0.45 µm	1 = 1/8" Hose Barb 2 = 1/4" Hose Barb 3 = 1/2" Hose Barb 4 = Luer Lock 5 = 1/2" Sanitary* 6 = 1/4" MNPT	1 = 1/8" Hose Barb 2 = 1/4" Hose Barb 3 = 1/2" Hose Barb 4 = Luer Lock 6 = 1/4" MNPT 7 = Side Bleed Valve	1 = 1/8" Hose Barb 2 = 1/4" Hose Barb 3 = 1/2" Hose Barb 3B = 1/2" Hose Barb with Filling Bell 4 = Luer Lock 5 = 1/2" Sanitary* 6 = 1/4" MNPT	1 = 1/8" Hose Barb 2 = 1/4" Hose Barb 3 = 1/2" Hose Barb 4 = Luer Lock 6 = 1/4" MNPT 7 = Side Bleed Valve	B = Buna E = EP S = Silicone V = Viton (or FKM) K = FFKM

Pre-Sterilized or Not
S = Pre-Sterilized
N = Not Sterilized

*When choosing the Sanitary Inlet/Outlet, the Luer Lock or bleed valve option is required for the Vent Port
** O-Ring is only available on Bleed Valve



One Chestnut Street
Nashua, NH 03060
603.880.4420
FAX: 603.880.4536

CriticalProcess.com

The information contained herein is subject to change without notice. The Critical Process Filtration logo is a trademark of Critical Process Filtration, Inc. Viton is a trademark of DuPont Performance Elastomers L.L.C.
© 2025 Critical Process Filtration, Inc. • All Rights Reserved

Data Sheet DPPC Micro DS Rev -