



PTR Micro Capsule Filters

PTFE Membrane

PTR Micro Capsule filters consist of a single layer of validated Polytetrafluoroethylene (PTFE) membrane and are used for sterilizing process gas applications and tank vents. The PTR Micro Capsule's membrane has a pore size of 0.22 μm .

PTR Micro Capsule filters are hydrophobic and have high air flow and low pressure drops. Each Micro capsule is individually integrity tested using the bubble point method before it is released from manufacture.

Critical Process provides unrivaled delivery times, technical consulting before purchasing, and very competitively priced high-performance products. Our comprehensive testing & analysis and validation services support your team whenever they need it. Your process experts partnering with our filtration experts is how we deliver your company's solution right the first time.

Sterilizing Filters

Tank Vent & Process Gas




MICRO CAPSULES – Nominal Dimensions

Body Length: 1.9 in. (4.8 cm)

Overall Length: 2.8 to 3.8 in. (7.1 to 9.7 cm)

Outside Diameter: 2.6 in. (6.6 cm)

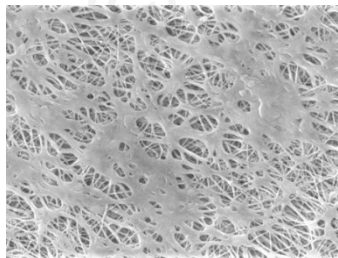

PTR Micro Capsule filters are recommended for:

- Tank Vents
- Compressed Air
- Pressurized Gases
- Fermentation Air

Maximum Operating Parameters

MICRO CAPSULES	
Gases Operational Pressure	60 psi at 68 °F (4.14 bar at 20 °C)
Operating Temperature	110 °F at 30 psid (43 °C at 2.07 bard)
Forward Differential Pressure	Gas - 60 psi at 68 °F (4.14 bar at 20 °C)
Reverse Differential Pressure	50 psid at 68 °F (3.45 bard at 20 °C)
Recommended Changeout Pressure	35 psid (2.41 bard)

Sanitization & Sterilization

Autoclave	250 °F (121 °C), 30 min, 5+ cycles
Chemical Sanitization	Performed using industry standard concentrations of hydrogen peroxide, peracetic acid, sodium hypochlorite and other selected chemicals.

Integrity Testing

PORE SIZE	BUBBLE POINT MINIMUM*	
	PSIG	BARG
μm		
0.22	18	1.24

* Bubble Point for membrane wetted with 60% IPA / 40% water

Filtration Area (Nominal)

Area	0.59 ft ²
	548 cm ²

Construction Materials

Filtration Media	Polytetrafluoroethylene (PTFE) Membrane
Media Support	Polypropylene
End Caps, Center Core, Outer Support Cage, Micro Capsule Housing	Polypropylene
Sealing Method	Thermal Bonding

Validation

PTR Micro Capsules are validated using test procedures that comply with the intent of BFE protocols for the determination of bacterial retention in filters used for air and gas filtration. The challenge level for the 0.22 µm filter membrane is 7.5×10^7 organisms per cm² of filter media of *Brevundimonas diminuta* (ATCC 19146).

PTR Micro Capsule filters are also validated using test procedures that comply with ASTM F 838-15(ae1) protocols for the determination of bacterial retention in filters used for liquid filtration. The challenge level is a minimum of 10^7 *Brevundimonas diminuta* organisms per cm² of filter media. CPF filters have > 7-log removal when challenged with the *Brevundimonas diminuta* (0.22µm meets the FDA definition of sterilizing grade filters).

Validation Guides available upon request.

Endotoxins

The levels of bacterial endotoxins in aqueous extracts from PTR Micro Capsule filters are below current USP limits as specified for water for injection.

Extractables

PTR Micro Capsule filters typically exhibit low levels of non-volatile residues.

TOC and Conductivity

The PTR Micro Capsule filters conform with TOC standards of USP <643> and the water conductivity standards of USP <645> after an appropriate flush with purified water.

Toxicity Compliance

Materials used to construct the PTR Micro Capsule filters are non-toxic and meet the requirements for the MEM Elution Cytotoxicity Test and the requirements for Biological Reactivity Tests in the current version of the United States Pharmacopeia (USP) for Class VI - 121 °C Plastics.


Non-Fiber Releasing

PTR Micro Capsule filters comply with Title 21 CFR sections 210.3 (b)(6) and 211.72, for non-fiber releasing filters.

FDA Compliance


Materials meet the requirements listed by the FDA as appropriate for use in articles intended for repeated food contact as specified in Title 21 CFR sections 174.5, 177.1500, 177.1520, 177.1630, 177.2440, and 177.2600 as applicable.

Water Flow Rates for
PTR Micro Capsules

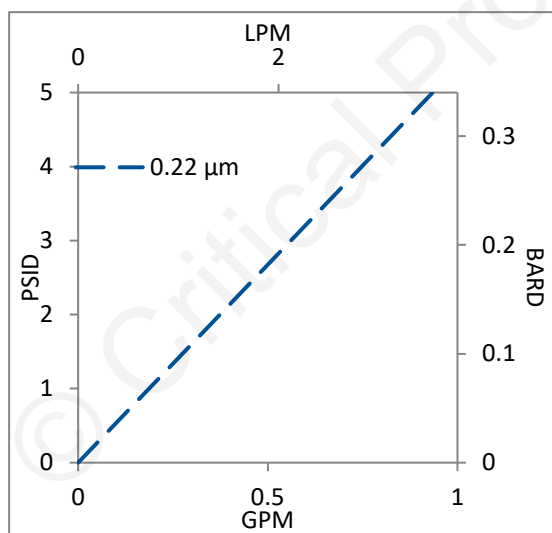

Flow rates for Micro Capsule filters are per filter. The test fluid is water at ambient temperature. Flows are tested using a Micro capsule filter with ½" Sanitary inlet and outlet ports. Rates will vary based on end configuration of the Micro capsule.

Air Flow Rates for
PTR Micro Capsules

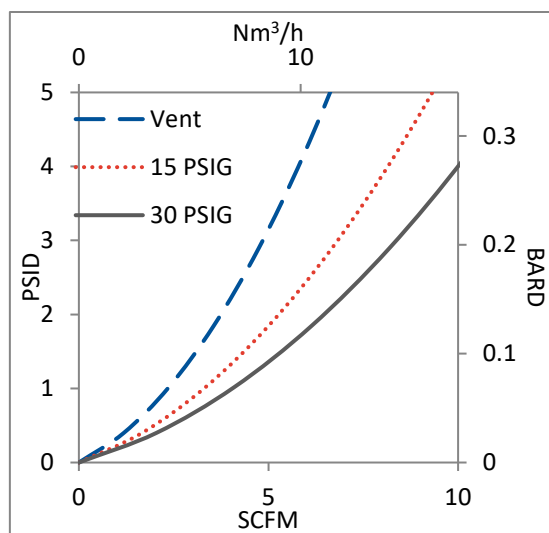

Flow rates for Micro Capsule filters are per filter. The test fluid is air at ambient temperature. Flows are tested using a Micro capsule filter with ½" Sanitary inlet and outlet ports. Rates will vary based on end configuration of the Micro capsule.

PTR Micro Capsule Filters Ordering Information

All Critical Process filters are configurable to meet customer specifications.
Fill in the corresponding codes in the boxes below to build your Part Number.

To consult with one of our technical team members, request a quote or place an order:
call (603) 880-4420 or [contact us here](#).

Micro Capsule Filters

MIC PTR N 000 -

Pore Size Code	Inlet	Inlet Vent Port	Outlet	Outlet Vent Port	Side Vent O-Ring**
-20 = 0.22 µm	1 = 1/8" Hose Barb 2 = 1/4" Hose Barb 3 = 1/2" Hose Barb 4 = Luer Lock 5 = 1/2" Sanitary* 6 = 1/4" MNPT	1 = 1/8" Hose Barb 2 = 1/4" Hose Barb 3 = 1/2" Hose Barb 4 = Luer Lock 6 = 1/4" MNPT 7 = Side Bleed Valve	1 = 1/8" Hose Barb 2 = 1/4" Hose Barb 3 = 1/2" Hose Barb 3B = 1/2" Hose Barb with Filling Bell 4 = Luer Lock 5 = 1/2" Sanitary* 6 = 1/4" MNPT	1 = 1/8" Hose Barb 2 = 1/4" Hose Barb 3 = 1/2" Hose Barb 4 = Luer Lock 6 = 1/4" MNPT 7 = Side Bleed Valve	B = Buna E = EP S = Silicone V = Viton (or FKM) K = FFKM

*When choosing the Sanitary Inlet/Outlet, the Luer Lock or bleed valve option is required for the Vent Port

** O-Ring is only available on Bleed Valve



One Chestnut Street
Nashua, NH 03060
603.880.4420
FAX: 603.880.4536

CriticalProcess.com

The information contained herein is subject to change without notice. The Critical Process Filtration logo is a trademark of Critical Process Filtration, Inc. Viton is a trademark of DuPont Performance Elastomers L.L.C.
© 2025 Critical Process Filtration, Inc. • All Rights Reserved

Data Sheet PTR Micro DS Rev -